


## **DBK - SEDIMENTATION SHAKER**

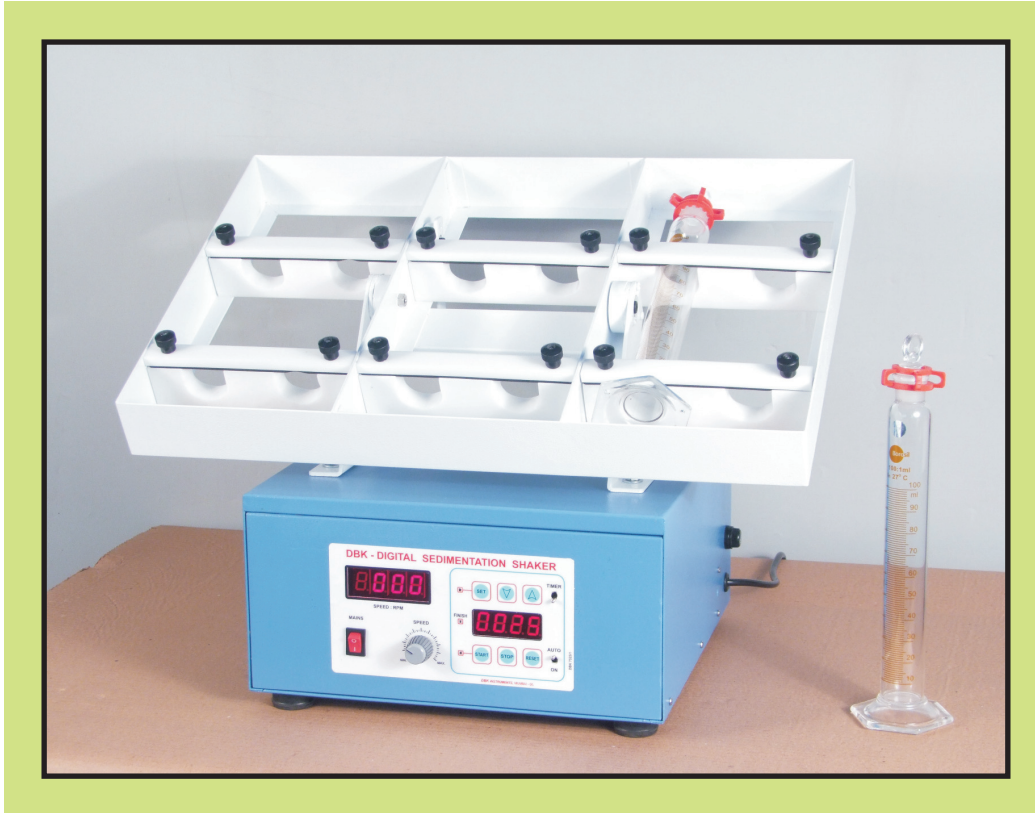

## **DBK - SEDIMENTATION SHAKER**

### **APPLICATION:**

Sedimentation shaker is device to determine the sedimentation value of the flour as a measure of protein quality & quantity, based on swelling of protein globules and agglomeration of starch into protein globules in the lactic acid medium.

### **FEATURES :**

Smooth movement at an angel of  $\pm 30^\circ$  from horizontal.  
Holding tray suitable for 50ml / 100ml cylinder.  
Easy mechanism for Loading & Unloading of cylinder

### **SPECIFICATIONS :**

Speed : 30 RPM (Fixed Speed / Digital Variable Speed).  
Cylinder : 50 / 100ml X 6 Nos.  
Rocking Angle :  $\pm 30^\circ$  from Horizontal (Fixed).  
Power Supply : 230V.  $\pm 10\%$  AC. 50Hz.  
Dimension : 480mm W x 480mm D x 350mm H approx

OPTIONAL : Digital timer 1 Min to 99Hr : 59 Min to switch off shaking.

## ***DBK INSTRUMENTS***

MFRS. : SCIENTIFIC LABORATORY EQUIPMENT'S AND ELECTRONIC INSTRUMENTS  
28, Interlink Ind. Estate, Caves Road, Jogeshwari (E), Mumbai-400 060. Tel.: (+91-22) 2837 5586, 3295 1914 Fax : 2838 1914  
Email : dbklab@hotmail.com Website : www.dbkinstruments.in